

Vials by Autosampler Model

Find autosampler vials the easy way!

Simply select your autosampler make and model. then show you the vials your autosampler requires!












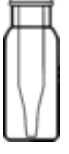
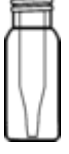






Agilent/HP
Alcott
Alltech
Beckman
Bruker
Carlo Erba/Fisons
Chrompack
Dani
Dionex/Gynkotek
Dynatech/Precision
Gerstel
Gilson
HTA
Jasco
Kontron
LDC
Leap technologies
Merck/Hitachi
Perkin Elmer
Pharmacia
Shimadzu
Spark Holland
Spectra-Physics/TSP
Tekmar-Dohrmann
Thermo Quest
Tosoh
Unicam
Varian
Waters




Vials by Pictures

12x32mm Autosampler Vials

Click on the picture of the vial you want...

	<u>Crimp Cap Vials</u>	<u>Snap/Crimp Vials</u>	<u>8mm Screw Cap Vials</u>	<u>9mm Screw Cap Vials</u>	<u>10mm Screw Cap Vials</u>
<u>Standard '2ml' Vials</u>					
<u>Champagne Vials</u>					
<u>X-Vials</u>			N/A		
<u>Other Micro-sampling Vials</u>					

Other Sizes

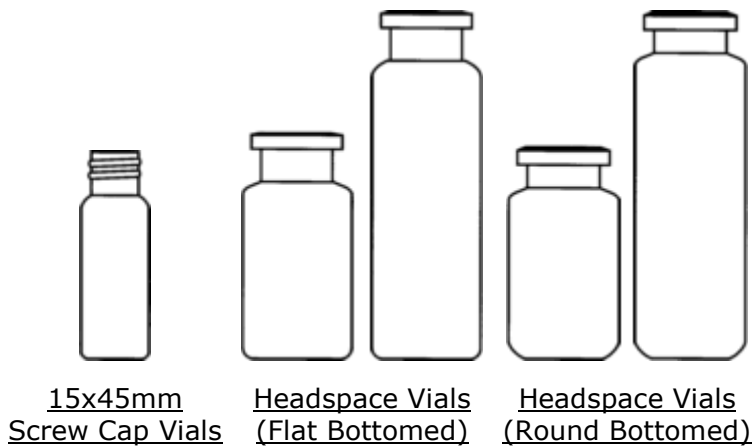

Vials by Closure:

Crimp Cap Vials



These vials are closed with 11mm aluminium crimp caps. Crimp caps provide the tightest seal, reducing the chance of sample evaporation. They are economical and easy to use with our vial crimping tools.

Snap/Crimp Cap Vials



Snap caps are often more convenient than crimps as they can easily be applied and removed by hand. Snap/crimp vials can be closed with either crimp caps or snap caps.

8mm Screw Cap Vials (Standard)



8mm (8-425) screw thread vials are the original, standard opening autosampler vials. Use with 8mm caps and septa.

9mm Screw Cap Vials



9mm screw top vials have a 9-425 GPI thread and can be used in all autosamplers that use 11mm crimp cap vials. Use with 9mm caps and septa.

10mm Screw Cap Vials (Wide Mouth)



10mm (10-425) screw cap vials feature a wider opening than other screw top vials. This increased target area makes sample preparation easier. It also reduces the chances of bent or broken needles during sampling. Use with 10mm caps and septa.

Vials by Size

12x32mm Vials



12x32mm vials are the most common type of autosampler vial and can hold 100 μ l-2ml of sample. They are available in the widest variety of materials and closure types.

15x45mm Vials



15x45mm vials can hold up to 4ml of sample. They are used by the Waters 48-position autosampler (and others).

Vials by Volume

Limited Volume Vials



Microsampling vials are standard 12x32mm vials with volumes less than 1.5ml. Be sure to look at our main X-vials page too!

1.5 - 2ml Autosampler Vials



1.5ml - 2ml 12x32mm vials are the most common variety of autosampler vial. They are available with the full range of closures.

>2ml Auto sampler Vials



Larger volume vials including headspace vials, sample storage vials and VOA vials.

International Co. for Lab. Instruments

45 Shalhoub st. Ain Shams, Cairo, Egypt Post Code 11311



+202 24916139 or +202 24916139



+202 24953441